


# Plasti-med®

"supporting healthy life"



**FDA** **CE**  
Registered Facility 1984

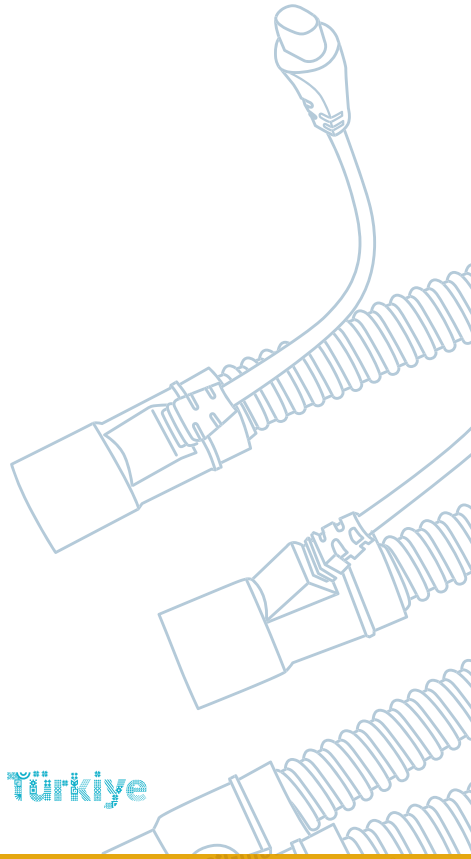

**Türkiye**

**HEATED CIRCUITS**




The spiral-heated wire provides humidified air to enhance comfort and reduce issues caused by dry air. It is designed to minimize condensation within the closed circuit.

Features that ensure ease of use include;

- Aspiration-friendly Y-connector,
- Siliconized temperature probe port,
- Temperature probe cable stabilizer and
- Adjustable humidification tube.

Manufactured from biocompatible materials in compliance with international standards.

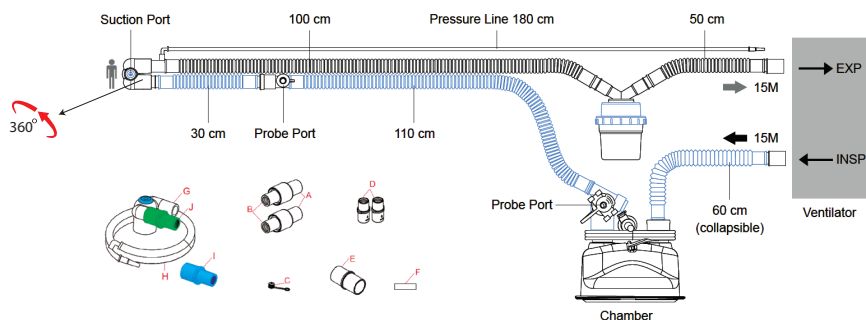


	Code	Y Connector	Elbow	Inhale	Exhale	Auto-Feed Chamber	Water Trap	Pressure Line
Infant	489071S	10 mm	15F with suction port	110+30 cm	150 cm	60 cm	✓	180 cm
	489061S	10 mm	15F with suction port	110+30 cm	140 cm	60 cm	-	180 cm
Adult	489011S	22 mm	22/15, 22 mm	150 cm	150 cm	60 cm	✓	-
	489001S	22 mm	22/15, 22 mm	150 cm	150 cm	60 cm	-	-



Plasti-med's new active humidification system, Embedded (spiral) Heated Wire Circuits, reduces condensation within the closed circuit and ensures the patient receives optimally heated and humidified air.

The integrated heated wires along the tubing eliminate the cold interface between the homogeneously heated respiratory gas and room temperature, preventing fluid buildup in the respiratory circuit.

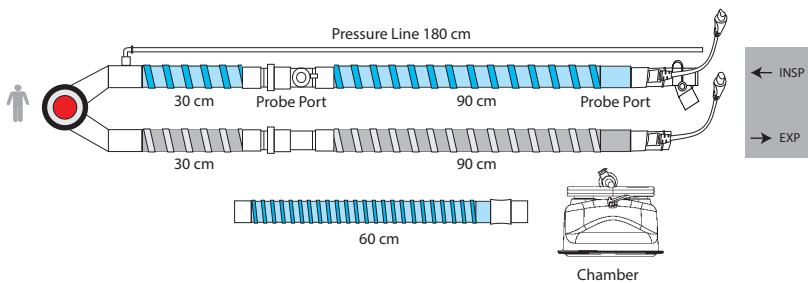

	Code	Y Connector	Elbow	Inhale	Exhale	Auto-Feed Chamber	Water Trap	Pressure Line
Infant	489071	12 mm	15F	90+30 cm	150 cm	60 cm	✓	180 cm
	489061	12 mm	15F	90+30 cm	120 cm	60 cm	-	180 cm
Adult	489011	22 mm	22/15, 22 mm	160 cm	160 cm	60 cm	✓	-
	489001	22 mm	22/15, 22 mm	160 cm	160 cm	60 cm	-	-

**EXPORT**  
ALL OVER THE  
WORLD

MEDICAL DISPOSABLES  
SINCE 1995 MANUFACTURER



Since 1995

## PLASTİ-MED PLASTİK MEDİKAL ÜRÜNLER SANAYİ VE TİCARET LİMİTED ŞİRKETİ

Deri Osb Mahallesi Yan Sanayi Cad. No: 13 Tuzla/İstanbul, TÜRKİYE

Tel: +90 (216) 591 04 91 Fax: +90 (216) 591 04 90

info@plasti-med.com

www.plasti-med.com

